





## TEST REPORT

(Electronic version)

ISSUE DATE: 2020-04-08

1J1号 F载报告



VERIFICATION WEBSITE: www.gttc.net.cn VERIFICATION CODE: ATOX-9073-34



No: 200068394 APPLICANT: JIAXING ZHENG QUN MEDICAL DEVICES CO., LTD.

BUILDING 3, XINJIA PARK, JIANGNAN RD. NANHU DIST. JIAXING CITY, ZHEJIANG ADDRESS:

SAMPLE INFORMATION:

THIRTY MEDICAL SURGICAL FACE MASKS

STYLE NO.: MF-ZQ2020.4

COLOUR: BLUE SIZE: 175×95mm

DATE RECEIVED/DATE TEST STARTED: 2020-04-06		
CONCLUSION:		
BACTERIAL FILTRATION EFFICIENCY	M	
APPEARANCE[3 PIECES]	M	
STRUCTURE AND SIZE[3 PIECES]	M	
NOSE CLIP[3 PIECES]	M	
MASK STRING[3 PIECES]	M	

NOTE: "M" –MEET THE STANDARD'S REQUIREMENT "F" –FAIL TO MEET THE STANDARD'S REQUIREMENT "---" –NO COMMENT

## REMARK:

THIS REPORT IS THE ENGLISH TRANSLATION VERSION OF THE REPORT 200068393.

ALL THE TESTED ITEMS ARE TESTED UNDER THE STANDARD CONDITION (EXCEPT FOR INDICATION).
COPIES OF THE REPORT ARE VALID ONLY RE-STAMPED.
THE EXPERIMENT WAS CARRIED OUT AT No. 1, ZHUJIANG ROAD, PANYU DISTRICT, GUANGZHOU, GUANGDONG, P.R.CHINA.

APPROVED BY:

ZiShan Guo SENIOR ENGINEER

郭子山



总部:广州市番禺区珠江路1号

花都实验室:广州市花都区狮岭镇旗岭河滨西路1号

电话:020-61994598/61994599 电话:020-37721161








## TEST REPORT

(Electronic version)

No:200068394

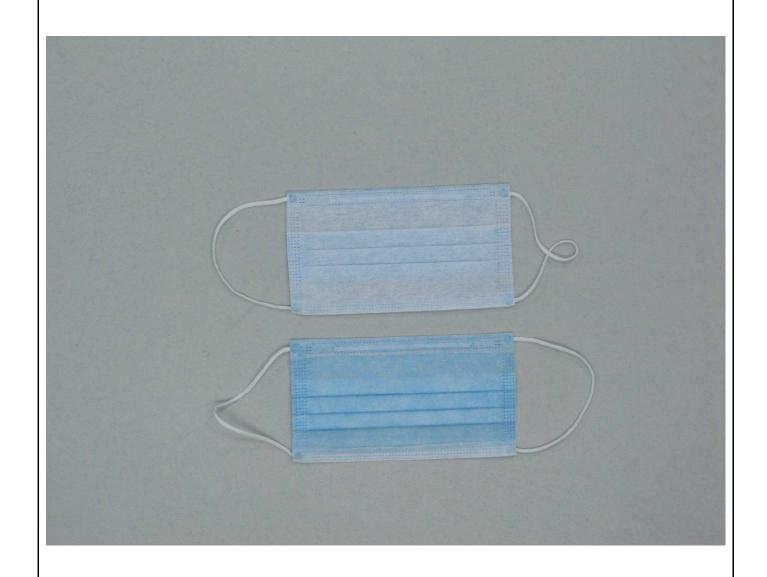



总部:广州市番禺区珠江路1号

花都实验室:广州市花都区狮岭镇旗岭河滨西路1号







## TEST REPORT

(Electronic version)

No: 200068394

BACTERIAL FILTRATION EFFICIENCY(%) 細菌ろ過効率(%)医療外科用マスク基準

(YY 0469-2011 ANNEX B, TEST BACTERIA: STAPHYLOCOCCUS AUREUS ATCC 6538, TEST AREA:  $40\text{cm}^2$ , FLOW RATE: 28.3L/min, MEAN PARTICLE SIZE: 3.0  $\mu$ m, RESULT OF THE POSITIVE CONTROL: 1.9 $\times10^3$  CFU, RESULT OF THE NEGATIVE CONTROL: <1CFU)

REQUIREMENT

 $\geq 95$ 

BFE<sub>1</sub> 98.5 BFE<sub>2</sub> 98.7 (YY/T 0969-2013)

BFE<sub>3</sub> 98.5 APPEARANCE[3 PIECES] 医療一般用マスク基準

(YY/T 0969-2013 5.1)

REQUIREMENT

**PASS** THE APPEARANCE OF THE MASK SHOULD BE

NEAT AND INTACT, AND THERE SHOULD BE NO

DAMAGE OR STAIN ON THE SURFACE.

(YY/T 0969-2013)

STRUCTURE AND SIZE[3 PIECES]

(YY/T 0969-2013 5.2)

REQUIREMENT

STRUCTURE: AFTER WEARING THE MASK, IT SHOULD BE ABLE TO COVER THE WEARER'S NOSE, MOUTH **PASS** TO JAW, AND SHOULD CONFORM TO THE DESIGN MAXIMUM SIZE DEVIATION(%):

SIZE, THE MAXIMUM TOLERANCE SHOULD NOT +1

EXCEED  $\pm 5\%$ . (YY/T 0969-2013)

REQUIREMENT

NOSE CLIP[3 PIECES]

(YY/T 0969-2013 5.3)

NOSE CLIP: THE PLASTIC NOSE CLIP SHOULD BE ATTACHED

TO THE MASK. THE NOSE CLIP LENGTH SHOULD PASS

NOT BE LESS THAN 8.0cm. NOSE CLIP LENGTH (cm): (YY/T 0969-2013)

1# 10.4 2# 10.4

3# 10.3

MASK STRING[3 PIECES]

(YY/T 0969-2013 5.4)

REQUIREMENT PASS

MASK STRING SHOULD BE EASY TO WEAR. THE BREAKING STRENGTH AT THE CONNECTION POINT BETWEEN EACH MASK BELT AND THE MASK BODY SHOULD NOT BE LESS THAN 10N.

(YY/T 0969-2013)



-End of Report----

总部:广州市番禺区珠江路1号

花都实验室:广州市花都区狮岭镇旗岭河滨西路1号

电话:020-61994598/61994599 电话:020-37721161